


tria

HAIR REMOVAL LASERS

Experience touchably smoother skin that uses diode laser technology preferred by dermatologists

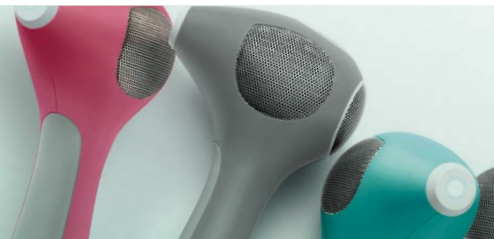

TRIA HAIR REMOVAL LASER TECHNOLOGY & BENEFITS


Hair Removal Lasers deliver permanent results for touchably smoother skin.

HOW DOES IT WORK

TRIA GETS TO THE ROOT OF HAIR REMOVAL

Tria Hair Removal Lasers are the only at-home hair removal devices featuring diode lasers, which disable the hair follicle by targeting the melanin at the root without damaging the surrounding skin.

1. TARGET THE HAIR FOLLICLE

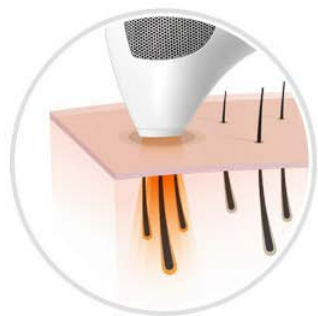


Tria lasers target the pigment in the hair follicle just like professional lasers.

2. FOLLICLE BECOMES DISABLED



The disabled hair begins to fall out, leaving behind more touchable, softer skin.


3. FOLLICLE PERMANENTLY DISABLED

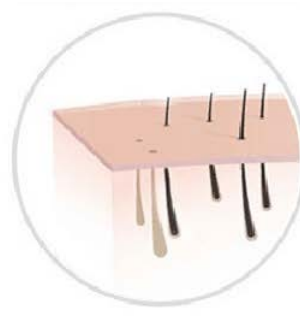

The laser permanently disables the hair follicles' ability to regrow hair.

TRIA VS. THE COMPETITION

Tria Hair Removal Laser 4X emits the highest level of energy of any at-home, FDA-cleared hair removal device.

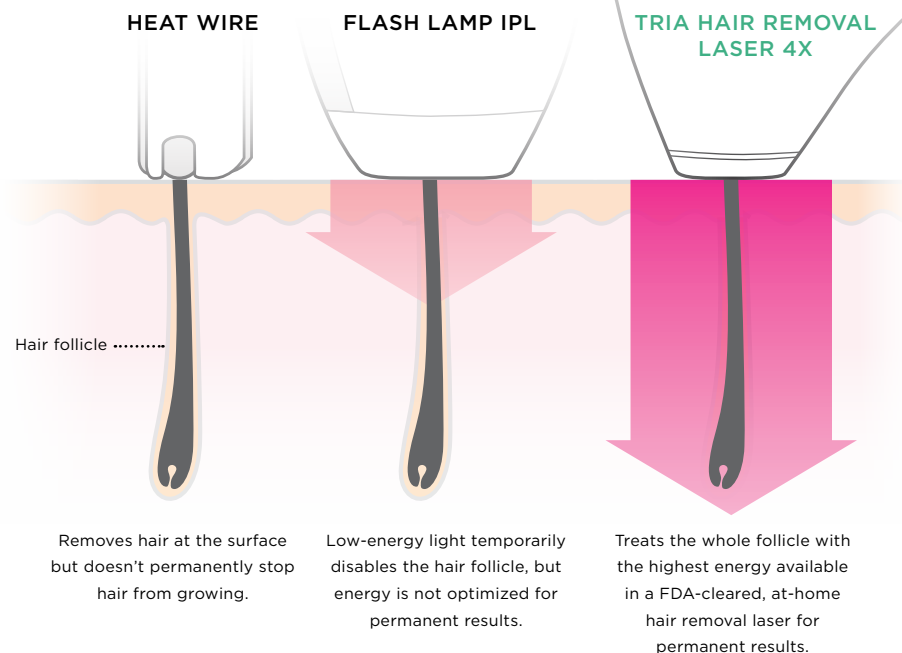

IS IT RIGHT FOR ME?

The Tria Hair Removal Laser 4X is designed to work on fair to medium skin tones and medium to dark hair—it will not work on gray, blonde or red hair.

SKIN TONE



White Ivory Beige Light Brown

Not intended for these skin tones—DO NOT USE



Medium Brown Dark Brown or Black

HAIR COLOR



Light Brown Medium Brown Dark Brown Black

Not intended for these hair colors—DO NOT USE



Blonde White or Grey Red

HOW TO USE

Three steps to touchably smoother skin.



STEP 1: PREP

Cleanse, shave and dry the area you want to treat. Prime your skin by applying Tria SmoothStart Calming Gel.

STEP 2: SET

Press the power button. Place skin sensor onto the skin you want to treat to test skin tone safety and unlock the device. Press the button to select your desired treatment level.

STEP 3: TREAT

Touch and hold the treatment tip flat against skin (correct contact signals light and a beep). Lift and place tip onto skin in an overlapping pattern.

